

a

PROJEKTARBEIT / STUDENT PROJECT

DUAL CHANNEL SINGLE MOLECULE LOCALISATION MICROSCOPY

concluded at the


TECHNISCHE UNIVERSITÄT (TU) WIEN

advised by

LUKAS VELAS

of

MORITZ SIEGEL

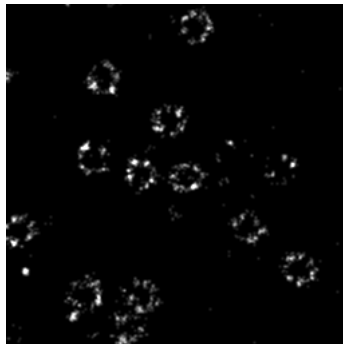

Vienna
May 25, 2021

moritz.siegel@tuwien.ac.at
<https://github.com/imrahilias/biophysics>

Contents

1	Introduction	2
2	Results	3
2.1	Phase Retrieval	3
2.2	Cramer Rao Lower Bound	9
3	Discussion	12
3.1	Phase Retrieval	12
3.2	Correction Collar	13
	List of Figures	14

Chapter 1

Introduction

Single Molecule Localisation Microscopy (SMLM) is a technique of fitting a full Point Spread Functions (PSF) to a stack of images containing reasonably spaced fluorescence signals.

To be able to accurately estimate the 3 dimensional location of the fluorescent molecule, the PSF must be known quite well. Simply put: If one knows the shape of a point source in varying degrees of defocus, one can guess the defocus and thus the z coordinate of the fluorescence molecule.

Chapter 2

Results

2.1 Phase Retrieval

It is possible to use different excitations to simultaneously measure different fluorescence markers, but for that the Point Spread Functions (PSF) has to be known for each wavelength. So the PSFs are estimated both for red and blue laser light, with the Phase Retrieval program by [Jesacher et al 21?]; each on a stack of 50000 STORM images of 100 nm beads stained with fluorescence dyes for both red (645 nm) and blue (488 nm) laser light.

To further complicate things, the new Olympus 1.5 NA objective comes with a Correction Collar to compensate for aberrations—which naturally vary slightly for both color channels. In order to find the best Correction Collar setting of the Olympus 1.5 NA objective for SMLM, the Point Spread Functions (PSF) is computed for each of the three settings of the correction collar {0.13, 0.17, 0.19}; using the program by [Jesacher et al 21?].

The Zernike modes of the PSF are shown for each of the three correction collar settings {0.13, 0.17, 0.19}, each for blue channel respective red channel in Figure 2.1 respective Figure 2.2.

For convenience the same results are additionally shown grouped by the three

correction collar settings {0.13, 0.17, 0.19}, now for both red and blue channel in Figure 2.3, Figure 2.4 and Figure 2.5.

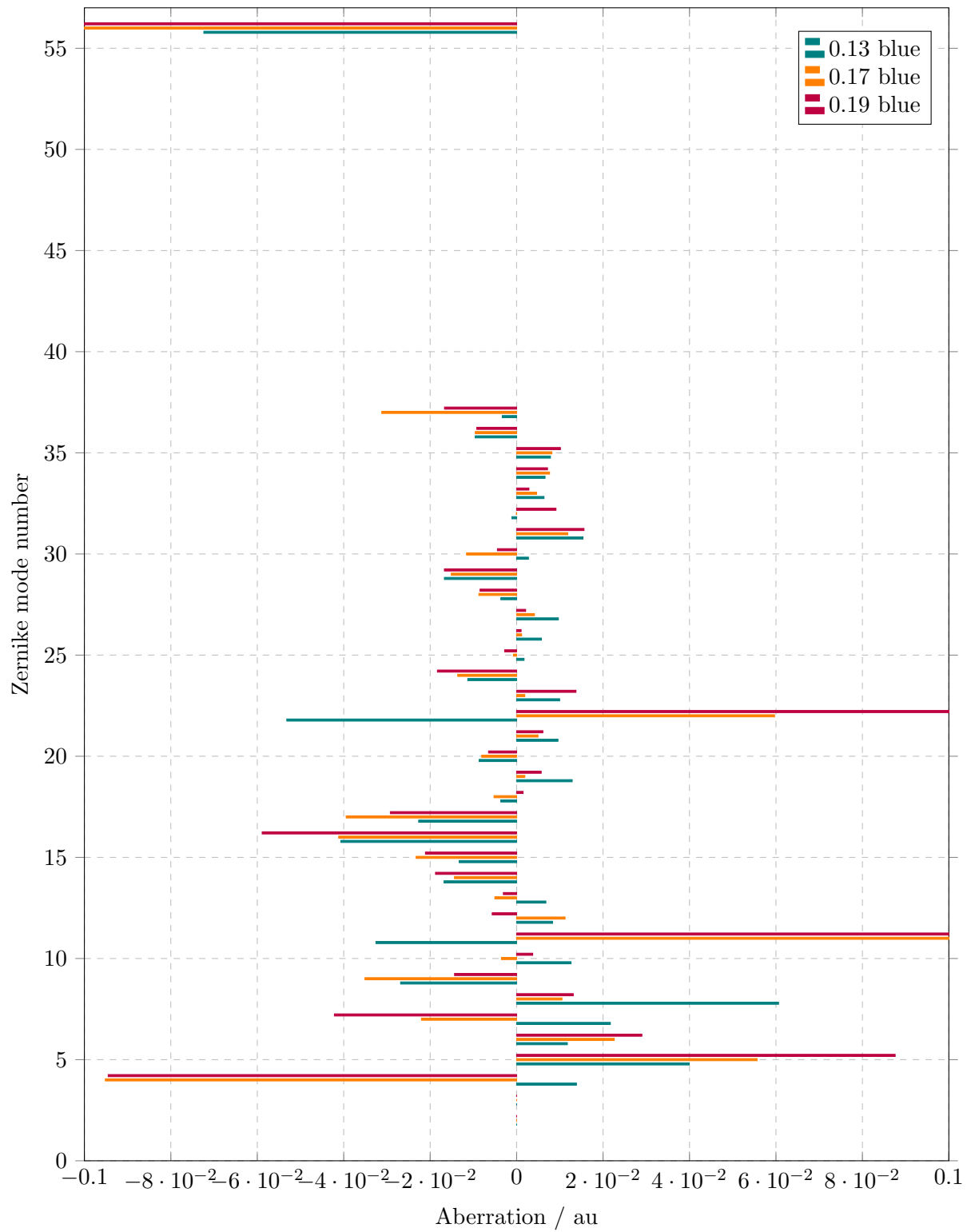


Figure 2.1: Blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for all three correction collar settings (0.13, 0.17, 0.19).

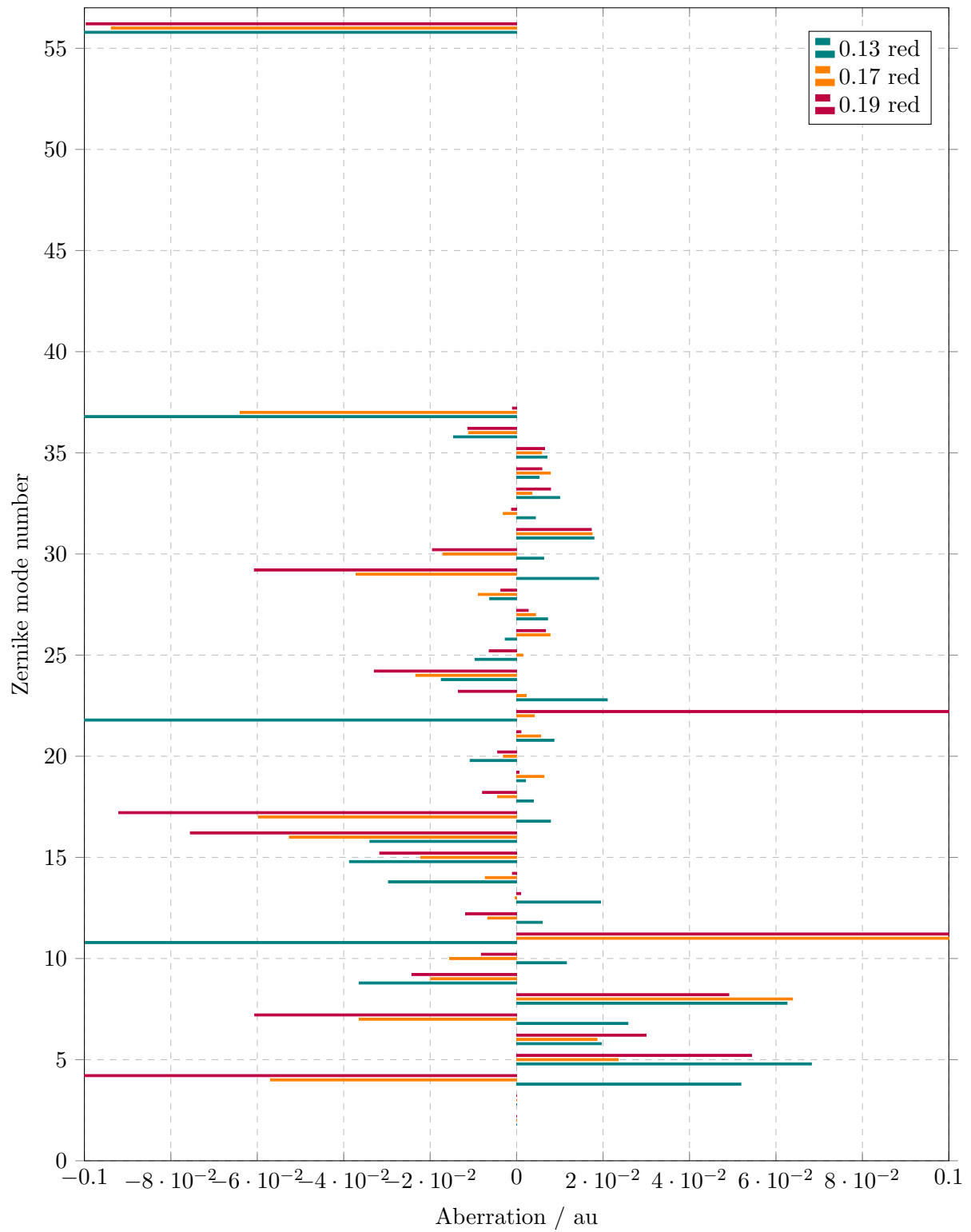


Figure 2.2: Red channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for all three correction collar settings (0.13, 0.17, 0.19).

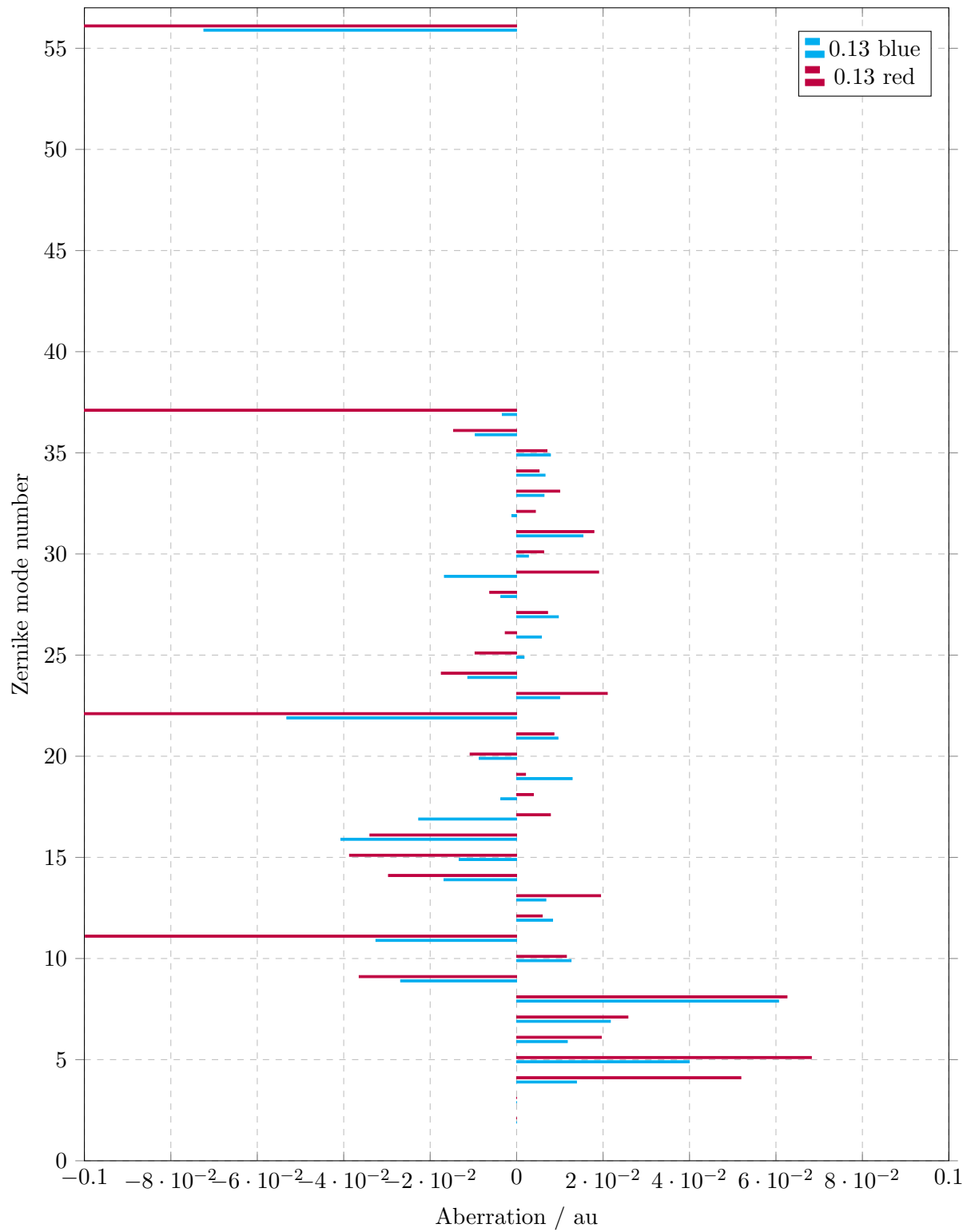


Figure 2.3: Red and blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for correction collar setting of 0.13.

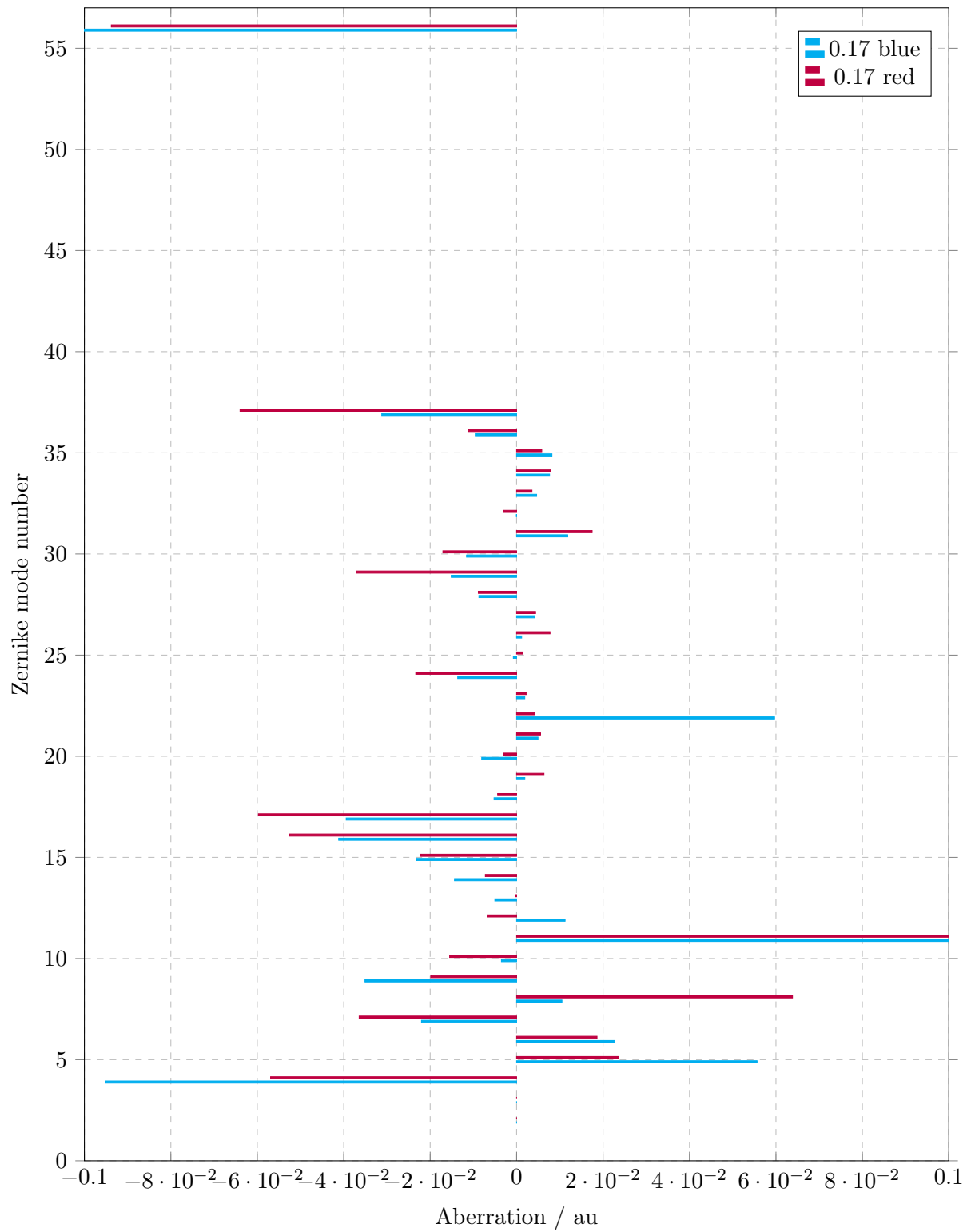


Figure 2.4: Red and blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for correction collar setting of 0.17.

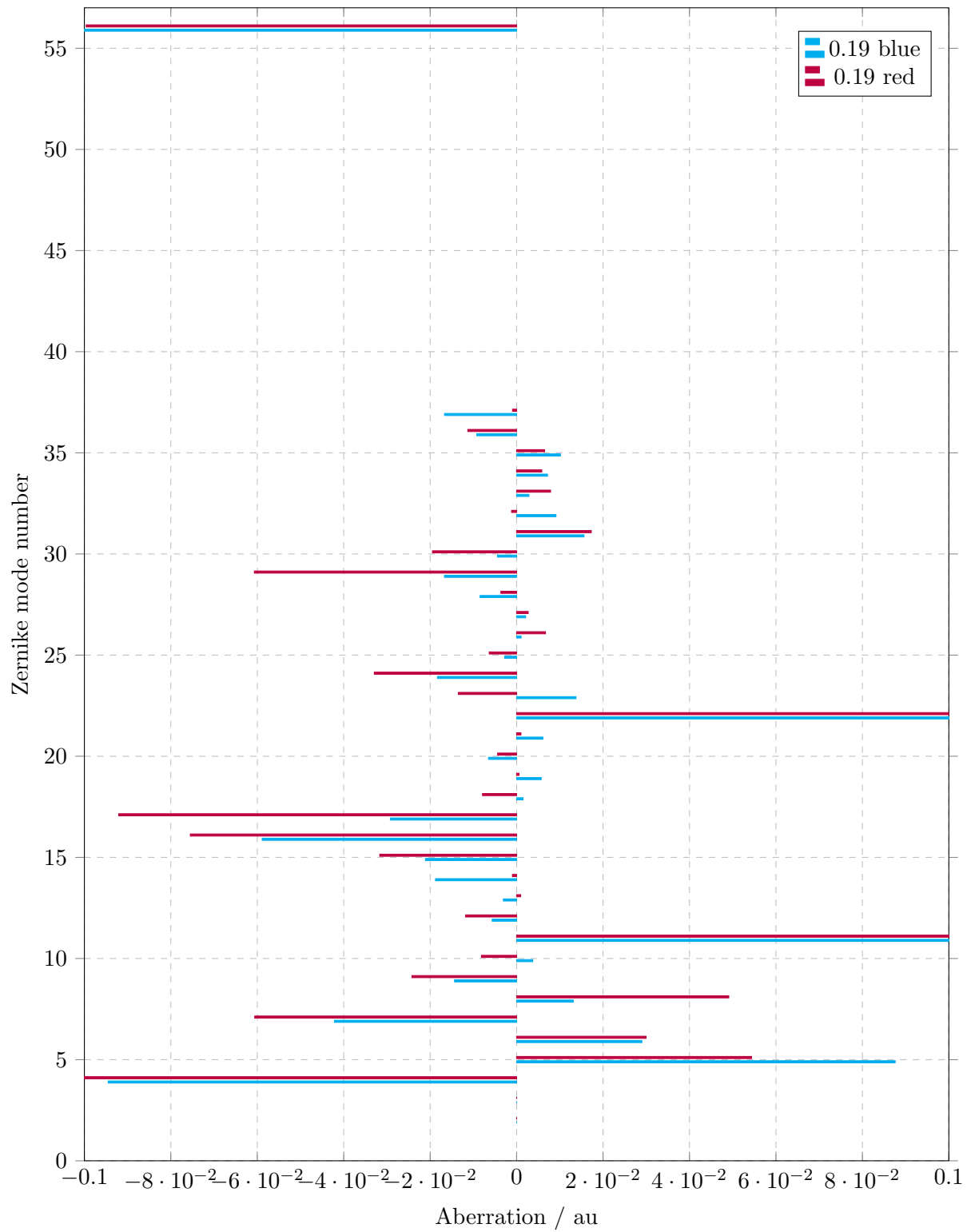


Figure 2.5: Red and blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for correction collar setting of 0.19.

2.2 Cramer Rao Lower Bound

In order to find the best Correction Collar setting of the Olympus 1.5 NA objective for SMLM, the Cramer Rao Lower Bound (CRLB) is computed using the program by [Jesacher et al 21?], with prior estimated Point Spread Functions (PSF) via phase retrieval [Jesacher et al 21?].

All CRLBs are computed at a defocus position of -500 nm, in steps of 5 nm from 0 to 500 nm. For all Estimations we assume identical arbitrary signal strength (2500) respective background (100).

The Estimated Cramer Rao Lower Bound for different correction Collar settings (varying linestyles for 0.13, 0.17, 0.19) of the Olympus 1.5 NA objective; for X, Y and Z axis (top, middle, and bottom); for both red and blue channel (colors) is shown in Figure 2.6.

Additionally plots grouped by X, Y and Z axis are shown in Figure 2.7, as well as grouped by different correction Collar settings in Figure 2.8.

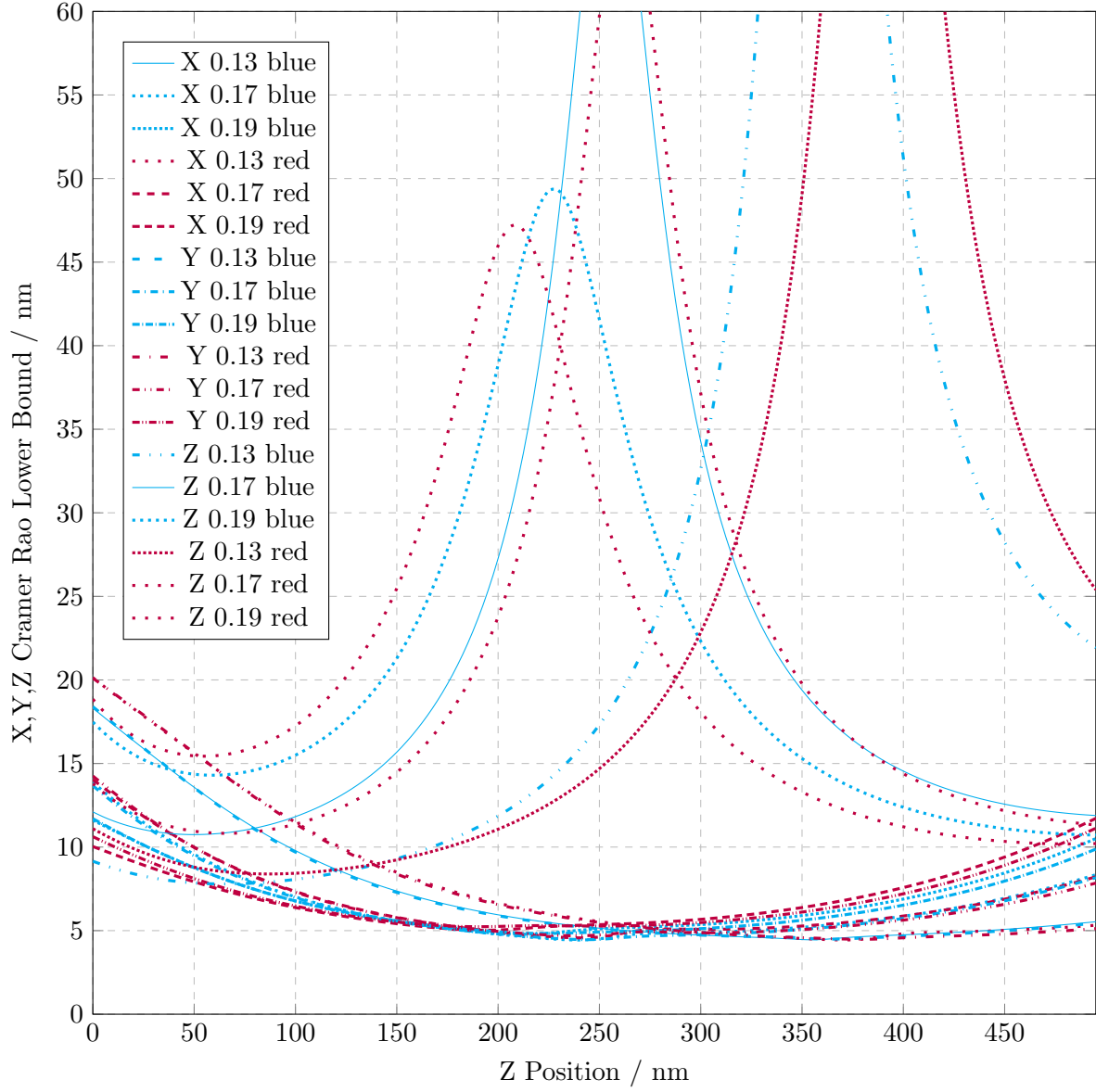


Figure 2.6: Estimated Cramer Rao Lower Bound for different correction Collar settings (varying linestyles for 0.13, 0.17, 0.19) of the Olympus 1.5 NA objective; for X, Y and Z axis (top, middle, and bottom); for both red and blue channel (colors).

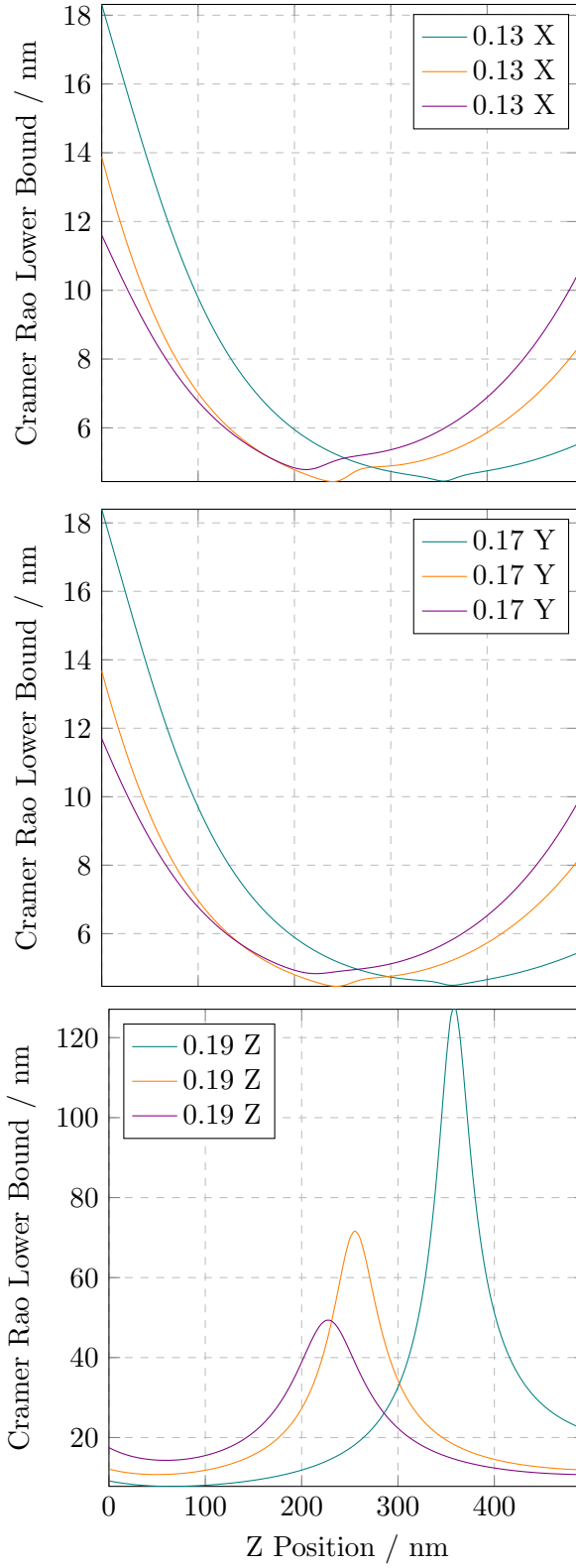


Figure 2.7: Estimated Cramer Rao Lower Bound for different correction Collar settings¹¹ (teal: 0.13, orange: 0.17, purple: 0.19) of the Olympus 1.5 NA objective; grouped by X, Y and Z axis (top, middle, and bottom).

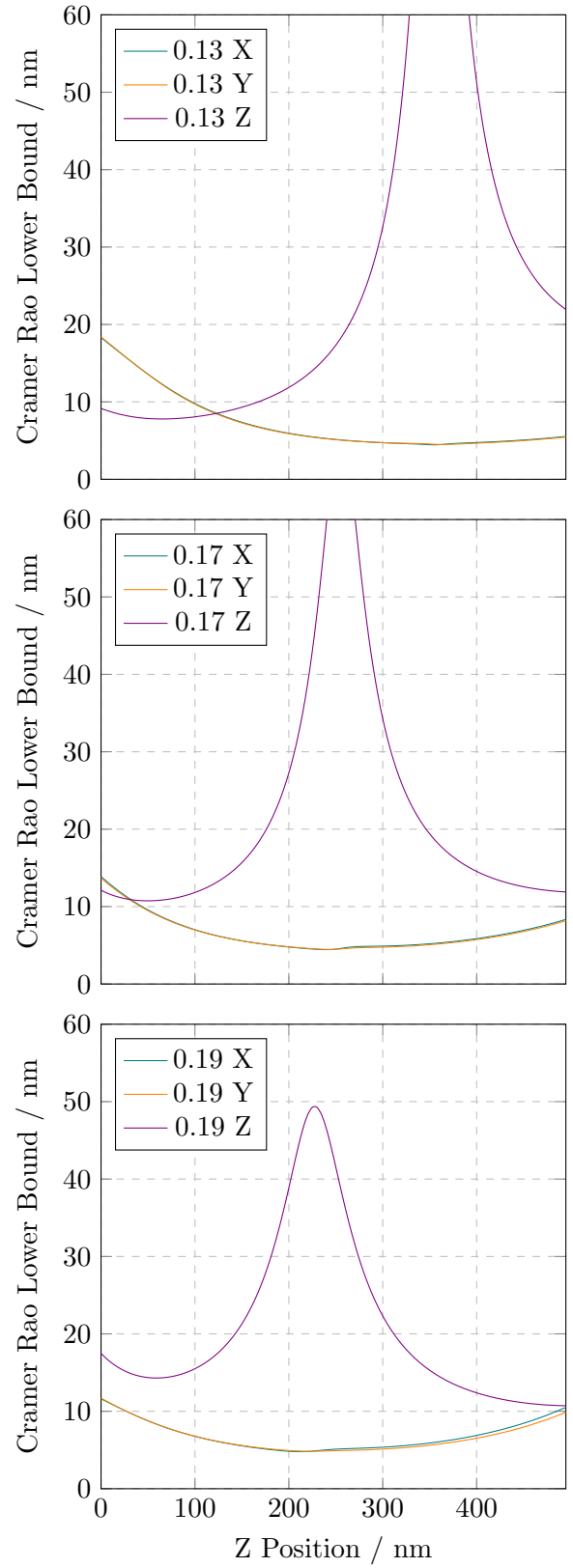


Figure 2.8: Estimated Cramer Rao Lower Bound for X, Y and Z axis (teal, orange and purple lines); grouped by different correction Collar settings (top: 0.13, middle: 0.17, bottom: 0.19) of the Olympus 1.5 NA objective.

Chapter 3

Discussion

3.1 Phase Retrieval

The image stacks look ok, even if not all are perfectly symmetrical—some errors are to be expected.

The estimated zernike modes of the PSF are mostly plausible; based on the modes and the order of magnitude.

This is further backed by the comparison of the results grouped by the three correction collar settings $\{0.13, 0.17, 0.19\}$, in Figure 2.3, Figure 2.4 and Figure 2.5. It is quite obvious, that the results for both red and blue channel are in the same order of magnitude—if not quite similar—for most of the Zernike modes, as suspected by theory!?

Yet the phase Retrieval program gave the error (high residual error for both the red and the blue channel when using the correction collar settings $\{0.17, 0.19\}$, is this a serious problem?

how can we avoid it?

3.2 Correction Collar

Based on the Figures 2.7 and Figure 2.8, and considering the fact—that we are interested in moderately defocusing (up to say 250 nm)—

one might conclude, that the most preferable setting for the correction collar is 0.13.

Yet the recommended setting for our microscope setup is 0.17! Which is backed by Lukas, based on the PSF stacks: One would have guessed the 0.17 is more sensitive to z since it changes much more with different z positions than 0.13.

Which is the one we should use?

List of Figures

2.1	Blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for all three correction collar settings (0.13, 0.17, 0.19).	4
2.2	Red channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for all three correction collar settings (0.13, 0.17, 0.19).	5
2.3	Red and blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for correction collar setting of 0.13.	6
2.4	Red and blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for correction collar setting of 0.17.	7
2.5	Red and blue channel Zernike modes $\{1 \dots 37, 56\}$ versus aberrations of PSF model via phase retrieval, for correction collar setting of 0.19.	8
2.6	Estimated Cramer Rao Lower Bound for different correction Collar settings (varying linestyle for 0.13, 0.17, 0.19) of the Olympus 1.5 NA objective; for X, Y and Z axis (top, middle, and bottom); for both red and blue channel (colors).	10
2.7	Estimated Cramer Rao Lower Bound for different correction Collar settings (teal: 0.13, orange: 0.17, purple: 0.19) of the Olympus 1.5 NA objective; grouped by X, Y and Z axis (top, middle, and bottom).	11
2.8	Estimated Cramer Rao Lower Bound for X, Y and Z axis (teal, orange and purple lines); grouped by different correction Collar settings (top: 0.13, middle: 0.17, bottom: 0.19) of the Olympus 1.5 NA objective.	11